


The PureBeam™ laser delivery system provides unprecedented flexibility, precision and visual control.

A variety of laser instruments have been available to ENT surgeons for decades. Today, after years of diligent research, development, and refinement, the PureBeam™ laser delivery system is enabling surgeons



to harness laser energy with unprecedented power, precision, and


results. The PureBeam™ system delivers a precision-focused laser spot size from Clinicon's extensive line of procedure-specific SureProbes™ with dissecting tips. The end result is greater flexibility and versatility, improved surgical access and visualization, and enhanced tactile control. Compared to current laser delivery systems, PureBeam™ is clearly the The Next Generation of CO₂ Laser Fiber Technology.

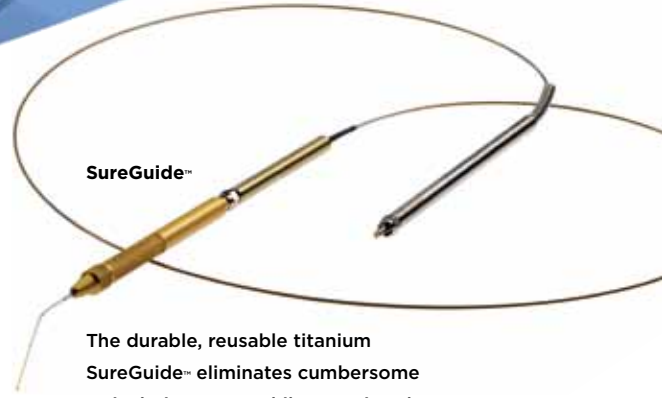

SureGuide™

The durable, reusable titanium SureGuide™ eliminates cumbersome articulating arms, while encasing the advanced, Clinicon fiber that preserves beam quality, intensity and efficiency.


The SureTouch™ universal handpiece is engineered to allow the quick replacement and interchangeability of the full line of Clinicon SureProbes™

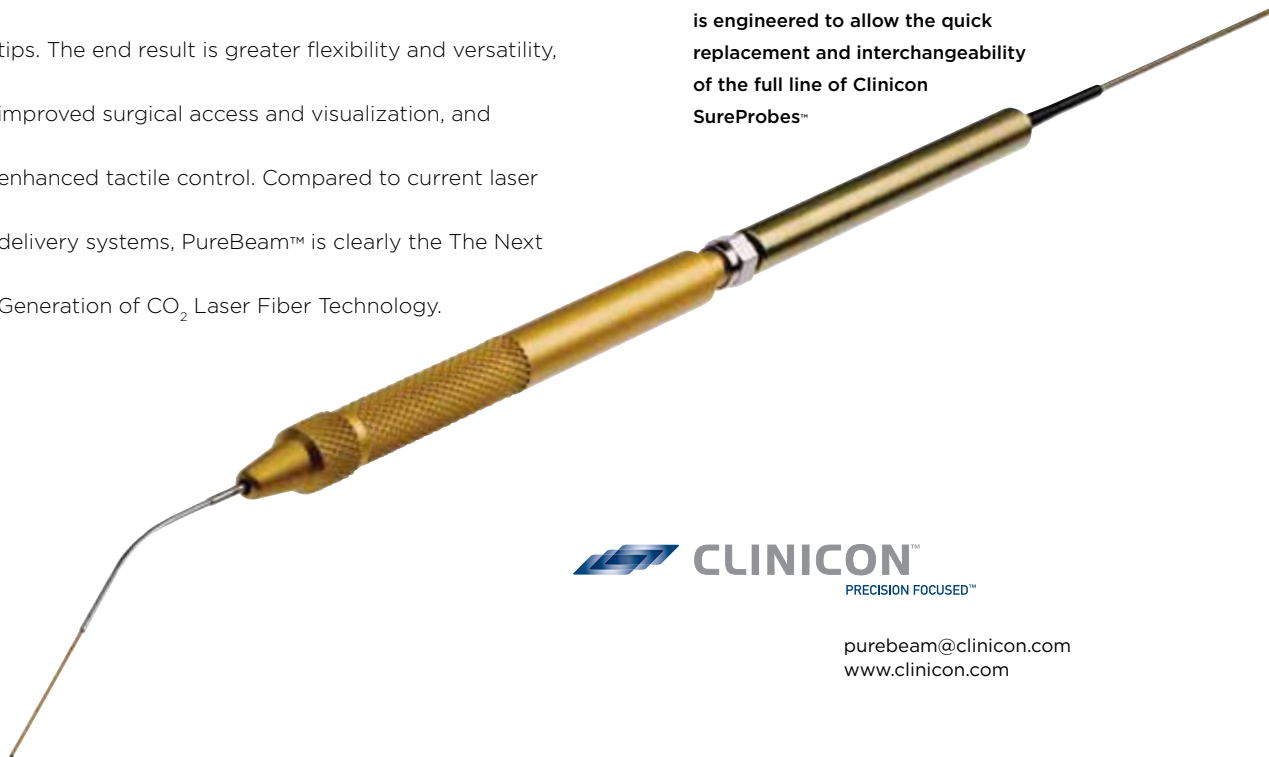

The Clinicon SureTouch universal handpiece and SureProbe tips are used to deliver CO₂ laser energy for the incision, excision, vaporization, ablation, coagulation or cauterization of soft tissue. Examples of some procedures in ENT are:

- Pulmonary Lesions
- Vocal Cord Dysplasia & Carcinoma
- Otosclerosis (Stapedectomy)
- Turbinectomy
- Bronchoectomy
- Cordectomy
- Polyps
- Larynx Leukoplakia
- Arytenoidectomy
- Papillomata
- Phenoidal Pituitary Tumor Removal
- Lymphangioma
- Hemiglossectomy
- Hypophysectomy
- Tuberoseptum Synechia
- Acoustic Neuroma
- Laser Assisted Myringotomy (LAM)
- Laser Assisted Uvulo-Palatoplasty (LAUP)
- In-excision Biopsy
- Removal of Benign & Malignant Lesions
- Frenectomy
- Gingivectomy
- Gingivoplasty
- Tuberosity Reduction
- Operculectomy
- Aphthous Ulcer Therapy
- Coagulation of Graft Donor Sites

SureProbes are available in a variety of straight and curved configurations.

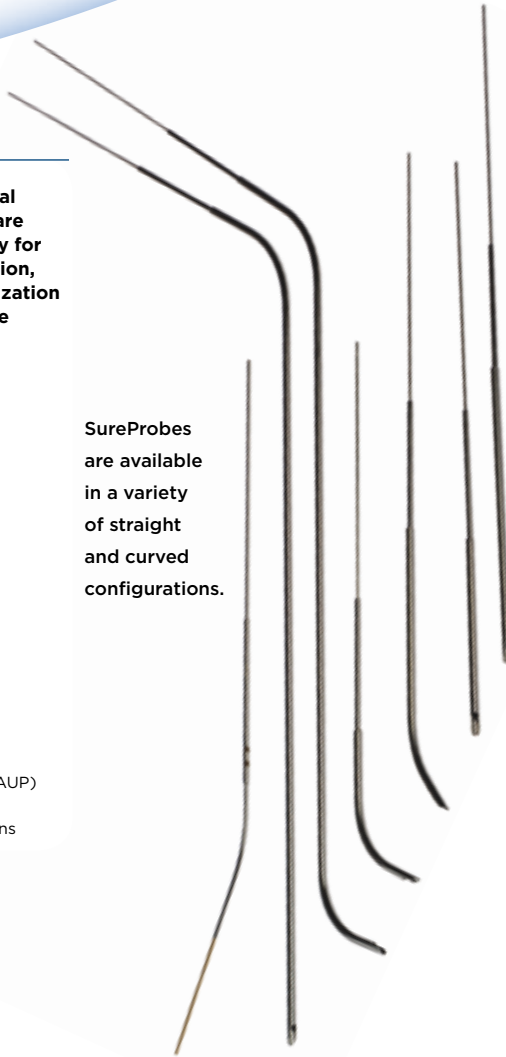

A scalloped dissecting tip with recessed fiber enables manipulation of tissue with the SureProbe without clogging the fiber aperture.



This allows the probes to be used in contact and non-contact mode while keeping them at an optimum distance from tissue for visualization during cutting and ablation. The probes are hollow, allowing gas flow to clear smoke from the surgical field.



All SureProbe ENT probes come packaged in a dispensing box that holds ten individually-sterilized probe pouches.

**FOR MORE INFORMATION OR TO PLACE AN ORDER CONTACT
1-800-CLINICON OR VISIT OUR WEBSITE AT WWW.CLINICON.COM**



Clinicon Corporation USA
3025-A Industry Street
Oceanside, CA 92054
T 800-CLINICON
T 760-439-1700
F 760-439-1798
purebeam@clinicon.com
www.clinicon.com

SUREPROBE™ ENT

PART # DESCRIPTION

4007	SureProbe™, ENT, 50mm Straight Dissecting Tip, Packaged Sterile Assembly
4003	SureProbe™, ENT, 50mm Single Curve Dissecting Tip, Packaged Sterile Assembly
4008	SureProbe™, ENT, 80mm Straight Dissecting Tip, Packaged Sterile Assembly
4004	SureProbe™, ENT, 80mm Single Curve Dissecting Tip, Packaged Sterile Assembly
4005	SureProbe™, ENT, 240mm Single Curve Dissecting Tip, Packaged Sterile Assembly
4006	SureProbe™, ENT, 240mm Double Curved Dissecting Tip, Packaged Sterile Assembly
4009	SureProbe™, ENT, BRONCH, 760mm Straight Probe, Packaged Sterile Assembly
4010	SureProbe™, ENT, BRONCH Ped, 380mm Straight Probe, Packaged Sterile Assembly
4011	SureProbe™, ENT, STAPES, 25mm Single Curve Probe, Packaged Sterile Assembly