


The PureBeam™ laser delivery system provides unprecedented flexibility, precision and visual control.

A variety of laser instruments have been available to surgeons for decades. Today, after years of diligent research, development, and refinement, the PureBeam™ laser delivery system is enabling surgeons to harness laser energy with



unprecedented power, precision, and results. The PureBeam™ system delivers a

precision-focused laser spot size from Clinicon's extensive line of procedure-specific SureProbes™ with dissecting tips. The end result is greater flexibility and versatility, improved surgical access and visualization, and enhanced tactile control. Compared to current laser delivery systems, PureBeam™ is clearly the The Next Generation of CO<sub>2</sub> Laser Fiber Technology.

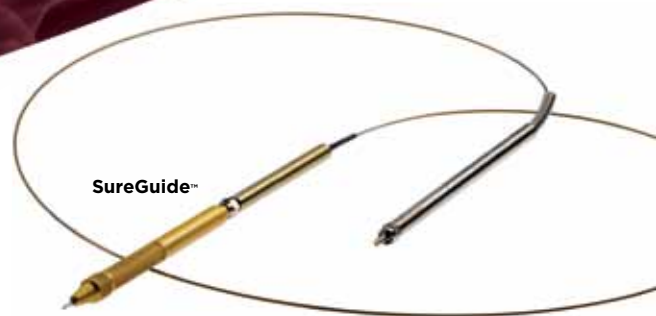

SureGuide™

The durable, reusable titanium SureGuide™ eliminates cumbersome articulating arms, while encasing the advanced, Clinicon fiber that preserves beam quality, intensity and efficiency.

The SureTouch™ universal handpiece is engineered to allow the quick replacement and interchangeability of the full line of Clinicon SureProbes™



**The Clinicon SureTouch universal handpiece and SureProbe tips are used to deliver CO<sub>2</sub> laser energy for the incision, excision, vaporization, ablation, coagulation or cauterization of soft tissue. Examples of some procedures in Spine are:**

Benign and Malignant Tumors  
 Astrocytomas  
 Lymphomas  
 Meningiomas  
 Chordomas  
 Sarcomas

SureProbes are available in a variety of straight and curved configurations.



All SureProbe Spine probes come packaged in a dispensing box that holds ten individually-sterilized probe pouches.



A scalloped dissecting tip with recessed fiber enables manipulation of tissue with the SureProbe without clogging the fiber aperture.



This allows the probes to be used in contact and non-contact mode while keeping them at an optimum distance from tissue for visualization during cutting and ablation. The probes are hollow, allowing gas flow to clear smoke from the surgical field.

## SUREPROBE™ SPINE

PART #	DESCRIPTION
4000	SureProbe™, Neuro, 40mm Straight Dissecting Tip, Packaged Sterile Assembly
4001	SureProbe™, Neuro, 50mm Straight Dissecting Tip, Packaged Sterile Assembly
4002	SureProbe™, Neuro, 80mm Straight Dissecting Tip, Packaged Sterile Assembly

**FOR MORE INFORMATION OR TO PLACE AN ORDER CONTACT  
 1-800-CLINICON OR VISIT OUR WEBSITE AT [WWW.CLINICON.COM](http://WWW.CLINICON.COM)**



**Clinicon Corporation USA**  
 3025-A Industry Street  
 Oceanside, CA 92054  
 T 800-CLINICON  
 T 760-439-1700  
 F 760-439-1798  
[purebeam@clinicon.com](mailto:purebeam@clinicon.com)  
[www.clinicon.com](http://www.clinicon.com)